

Cirrus™ Electric Ceiling Lift

iQ Technology™ and optional Integrated Scale provide smarter lifting, easier patient weighing

The new **Cirrus™ Electric Ceiling lift** combines two innovative features - an optional integrated patient weigh scale and iQ Technology - to create the smartest ceiling lift on the market. The integrated scale saves up to 12" in lifting height making it possible to accommodate facilities with low ceilings. A 450lbs, 600lbs and 750lbs model is available.

Tollos™ revolutionary iQ Technology™ is integrated into the lift to provide intuitive diagnostics for a smarter lifting program, including supplying data on number of lifts, hours in use, low battery operation, last/next service and service alerts every 1,000 lifts.

Our true six-point spreader bar, with optional locking carabiner clips, allows for maximum patient safety and can be used for effective in-bed repositioning - a task that older, two-point spreader bars just cannot effectively perform.

Plus, a new optional patent-pending wipedown lifting strap allows for the Cirrus™ to be completely disinfected, unlike older fabric lifting straps.

This lift is even further enhanced by our Curtain Track Pass Through option. This solution enables you to avoid cutting privacy curtains by allowing your existing privacy curtains to work seamlessly with the ceiling lift tracks.




The Cirrus™ Electric Ceiling Lift combines two innovative features - an integrated patient weigh scale and iQ Technology™ - to create the smartest ceiling lift on the market.



Optional wireless hand control - access all iQ Technology™ with convenient digital hand control.



Our true six-point spreader bar, with optional locking carabiner clips, allows for maximum patient safety and can be used for effective in-bed repositioning - a task that older, two-point spreader bars just cannot effectively perform.

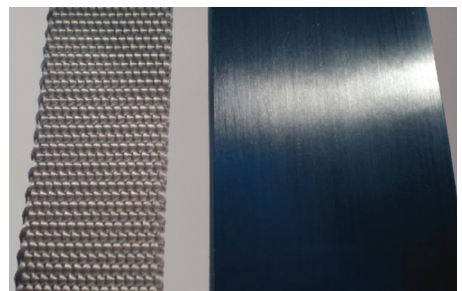

Tollos™ Infection Control lifting strap (right) can be wiped down to help prevent transmission of microbes vs. a traditional porous lift strap (left).

Cirrus™ Electric Ceiling Lift

Cirrus™ Electric Ceiling Lift

- Designed specifically for the institutional market, the new Cirrus™ Electric Ceiling Lift features a high output lifting motor and efficient gear drive for superior safety, durability and battery life.
- Available with optional continuous in track charge.
- Packed with an extensive array of safety and ease-of-use features.

iQ Technology™

Lifting and transferring has never been more intuitive with the integration of iQ Technology™. By utilizing proprietary electronics, mated to a digital readout and optional integrated scale, caregivers and maintenance staff can now track critical lift and service data by the push of a button. Data displayed includes total lifts, total hours, total battery charges, current battery level, total low battery warnings, last/next service plus visual alerts every 1,000 cycles for preventative maintenance.

Optional Integrated Scale

- Optional integrated scale saves up to 12 inches of lighting height and 40% cost versus external scale
- Eliminates maintenance including replacing batteries on external scale

Optional Hand Control with Digital Read-Out

- Access all iQ Technology™ with convenient digital hand control
- Access patient weight if using integrated scale
- Optional wireless hand control

Features

- 450 lbs, 600 lbs & 750 lbs lifting capacity
- iQ Technology™
- Integrated weigh scale (optional)
- Quick release trolley
- Remote hand control with display (optional)
- Wireless hand control (optional)
- True six-point spreader bar
- Optional patent-pending wipedown lifting strap
- Curtain Track Pass Through solution
- Controlled emergency lower without tools

TECHNICAL SPECIFICATIONS

Power Source:	2 x 12VSLA 2.7 Amp Batteries
Safe Working Load:	Model 450: 450 lbs (204 kgs) Model 600: 600 lbs (273 kgs) Model 750: 750 lbs (340 kgs)
Height (excluding trolley):	4.5" (11.4 cm)
Width:	12.5" (31.75 cm)
Length:	12.5" (31.75 cm)
Spreader Bar:	Six point
Spreader Bar Length:	23" (71.1 cm)
Spreader Bar Width:	15" (61cm)
Spreader Bar Depth:	10" (25.4cm)
Lifting Strap Length:	8' (243 cm)
Charger:	120VAC/24VDCCSA/UL Approved
Lifting Motor:	24VDC
Overall Weight:	27.2 lb (12.3 kg)
Charging:	End of Track or CC available
Emergency Stop:	Standard
Visual Battery Indicator:	Standard
Audio Low Battery Indicator:	Standard
Emergency Lowering:	Standard - Mechanical & Electric
Anti Free Wheel Safety System:	Standard
Up Limit Switch:	Standard
Down Limit Switch:	Standard
Remote Hand Controller:	Standard
Warranty:	1 Year
Integrated Weigh Scale:	Optional
Wishbone Style Spreader Bar - with Power Positioning	Optional
Compliant with:	ISO 13485 EN60601-1
CE	
UL/CSA	Motor, brake, charger & batteries

Formerly T.H.E. Medical

U.S. Corporate Office
1 Easter Court, Suite J
Owings Mills, MD 21117
410.363.1515 PHONE
888.220.4566 TOLL FREE

Canadian Office
75 Dymont Road
Barrie Ontario L4N 3H6
705.733.0022 PHONE
888.363.7224 TOLL FREE