

PERFECT FOR PATIENTS
PERFECT FOR CAREGIVERS




CREATING INNOVATION AND DIFFERENCE

PulmoAir

PulmoAir combines the benefits of lateral rotation, Low Air Loss and Pulsation into a single system. The system provides two selectable turning angles to bilaterally rotate the patient to a full or partial turn to reduce the risk of complications associated with immobilization. PulmoAir is also designed with microprocessor for simplicity, reliability and ease of use.

PulmoAir Mattress System is suitable for medium to high risk patients who are at risk of developing Pressure Ulcers stage 3 and 4.

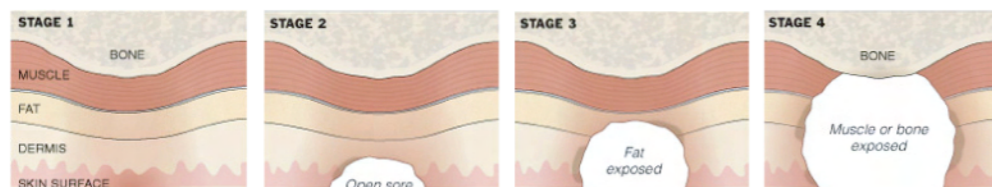

Various top covers are designed and available to satisfy your individual requirement in terms of therapy, prevention and patient comfort.

Mattress Control Unit

- Dimension: (L)12.3"x (W)6.7"x (H)4.9"
- Weight: 14.9 lbs
- Rated Voltage: 110-120 volts
- Max Current: 5A
- Rated Frequency: 60Hz
- Cycle Time: 3~95 mins
- Min/Max Pressure: 15~40 mmHg +/- 5 mmHg
- Static / Pulsation Function
- Protection Type: Class I, Type BF
- 10 Digital Adjustable Comfort Control
- Bottom Out Alarm: Visual / Audio
- Leaving Bed Alarm: Visual / Audio
- Power Failure Alarm: Visual / Audio
- User-Friendly LCD Control Panel
- Fowler Function
- Auto Function
- Selectable Rotation Angles
- Pulsation

Mattress Replacement

- Model No: PulmoAir
- Dimension: (L)80"x (W)36"x (H)10"
- Weight: 39.6 lbs
- Weight Capacity: 600 lbs
- Cell Material: Blue Nylon with PU Coating
- Cell Number: 18 cells
- Cover Material: Blue Nylon with PU coating finish
- Base Material: Black Polyester with PVC Coating
- CPR Quick Release System

Toll Free:
877-367-5507
TAMPA, FL OFFICE
TEL: 813-854-5464



PulmoAir



PulmoAir

provides you, the carer, with peace of mind
and your patient, with world class care & comfort


Perfect for patients and caregivers




Low slip base cover is a plus to keep the mattress in place on the bed



Low Air Loss manages moisture and temperature

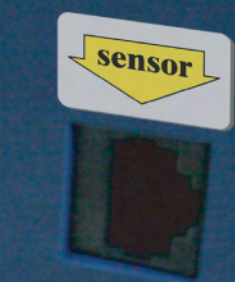

Incorporates with automatic sensing technology for bottoming out prevention and pressure adjustment

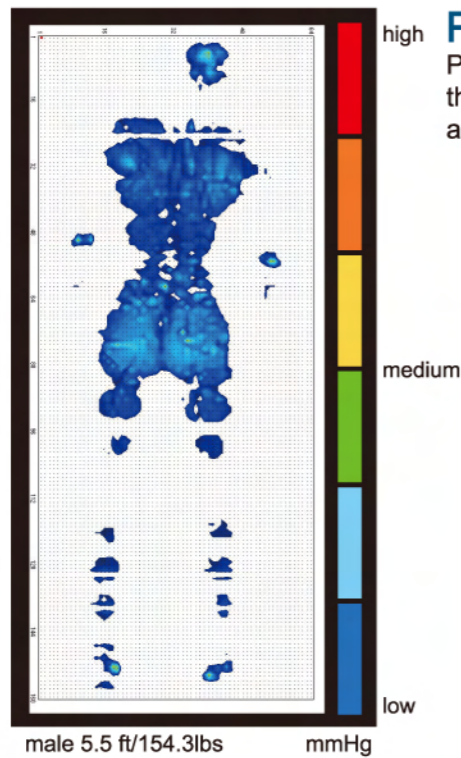

PRESSURE MAPPING

PulmoAir provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.



- The PulmoAir mattress replacement can accommodate patient's weight up to 600lbs
- Quilted top cover is highly vapor permeable and oversized providing a low friction and low shear against complications due to moisture
- Convoluted foam base provides for additional safety for patient
- Modular design ensures ease of use and cleaning
- Microprocessor controlled for simplicity, reliability and ease of use
- Provides optimum comfort by patient's weight and height
- Two selectable rotation angles with Low Air Loss

PulmoAir is your partner in protection and treatment of pressure ulcer



Cable management keeps electrical cables off the floor



CPR quick release system for rapid deflation

The Low-Pressure Alarm is triggered when the air pressure of the mattress is lower than the default setup.

Patient Care provides immediate firm surface during nursing procedures.

10 adjustable comfort settings provide different levels of comfort depending on the patient's therapy requirement.



The POWER FAILURE LED is triggered during power failure, providing a safer environment for patients.

Auto profile sets mattress pressure based on patient's height and weight.

Pulsation provides alternating pressure therapy.



Hanging bracket fits on all type of beds